

Alrick.

2300 Tilt Partner Bed Instruction Manual

Product Code: 7-23TP

Please Carefully read instructions before using the bed.



INTRODUCTION

This instruction manual contains information on the instructions for use and explanations of the signs and symbols used. These instructions for use contain instructions for operating the medical care bed.

These instructions for use may contain inaccuracies or printing errors. The information provided here will be updated periodically and changes due to product maintenance will be included in later editions. Changes or improvements are possible at any time without prior notice. If you have any questions, please contact our customer service team.

The instructions for use must be read and used by every person who operates the medical bed.

PURPOSE

The intended purpose of the Alrick 2300 Partner Bed is to be used for patients in hospital and healthcare facilities. These beds are electrically powered to allow the bed to be placed in various positions suitable to the patient. It is a medical device intended for human medical use only.

The Alrick 2300 Partner Bed moves on castors and features infinitely adjustable motorised height settings. The bedding area can be adjusted by multiple motors, allowing for positions such as sitting or exiting position.

The intended user of the Alrick 2300 Partner Bed is any bariatric (but not limited to) individual who has a disability or a health condition.

The bed should only be used by individuals trained or experienced in proper handling. Users must be instructed on its correct use and familiarize themselves with these instructions. Common users include nursing staff, physicians, therapists, family members, technicians, residents, and cleaning staff.

SYMBOLS EXPLAINED



Manufacturer



For Indoor Use (house)



Type B Application Part
(classification according to the degree of protection against electric shock)



Product with a Thermofuse



Safety Isolating Transformer




Device Protection Class II
(classification according to the type protection against electric shock)



Demko Approval

PRODUCT SERIAL NUMBER PLATE

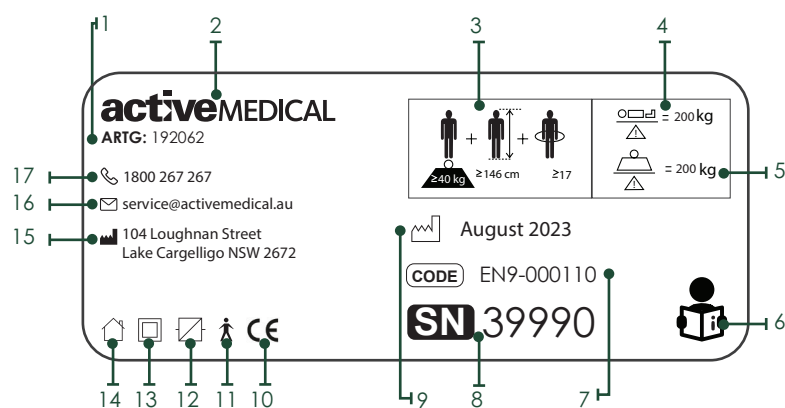

Serial Number Plates

SERIAL NUMBER PLATE EXPLAINED

- 1 ARTG Number
- 2 Company Logo
- 3 Adult Identification (minimum weight, minimum height, minimum BMI)
- 4 Maximum Occupant Weight
- 5 Maximum Safe Working Load
- 6 Read Instruction Manual - Follow Instructions For Use
- 7 Product Code
- 8 Product Serial Number
- 9 Manufacture Date
- 10 CE Marking
- 11 Type B Application Part
(classification according to the degree of protection against electric shock)
- 12 Switching Power Supply Unit
- 13 Device Protection Class II
(classification according to the type protection against electric shock)
- 14 For Indoor Use Only
- 15 Manufacturer Address
- 16 Company Service Email Address
- 17 Company Phone Number

PRE OPERATION CHECK

With the bed in position, check the following:

Electric beds should ONLY be used where earth leakage is fitted to the power source.

- Check power cord is plugged into 240V supply and is active.
- Check power cord is clear of castors and any moving part of bed.
- Ensure power cord has enough slack so cord is not strained when the bed is raised and lowered.
- Position bed at least 50mm from the wall in the lowest position. Lock castors. This position is critical for not damaging the wall when raising the bed.
- Before using bed, perform all functions to ensure they are working properly. If not please call your service technician or call the Active Medical Technician Team on 1800 282 001

** When fitting safety sides to this bed, please beware of the possible crushing points. The clamps can create crushing points on the rails. **

2300 PARTNER BED - MODELS AVAILABLE

2300 Partner King Double (900mm each side)

SPECIFICATIONS

Size	King Double
Height Range (mm)	360 - 800
Backrest Angle (degrees)	77
Kneebreak Angle (degrees)	40
Trendelenburg Angle (degrees)	16
Reverse Angle (degrees)	16
SWL (kg)	300 (each)
Overall External Dimensions (mm)	2050 x 1830
Sleeping Surface (mm)	2000 x 900 (each side)
Mattress Size (mm)	1960 x 870 (each)
TGA Number	ARTG 192062
Manufacturer 	Active Medical 104 Loughnan Street, Lake Cargelligo, NSW 2672

ELECTRICAL DATA

U Inv	230V ~, 50Hz
I In	Max. 1.5A
U Out	24V ... , max. 70W
Op.	10%, Max. 2min./18min.
Protection Class	IPx4
Insulation Class	II Type B
Noise	48DBA



For indoor use (house)



Safety isolating transformer



Type B equipment, as per EN60601+1



Device Protection Class II



Product with a thermo fuse




Demko approval

RE-Tuv Approval

IMPORTANT

- Check electric cords regularly and notify an electrician or Active Medical on 1800 282 001 immediately. If any damage occurs, cease to use and remove mains power until checked and repaired.
- No wiring is to be modified without authorisation from the manufacturer.
- Ensure the small extension cable is correctly connected to the bed socket and to the 240 volt supply when the bed is in a stationary position. A small light on the control box will illuminate showing a power supply.
- Low Low model beds should be left in the lowest position when a nurse is not in attendance.
- Nurses, cleaners or anyone using Alrick Electric Beds must make sure that the main power cord is clear of any moving parts on the bed before turning the power on.
- Electric Beds should ONLY be used where earth leakage is fitted to the power source.
- When patient condition could lead to patient entrapment, the mattress support platform should be left in the flat position when unattended.

OPERATOR CONTROLS

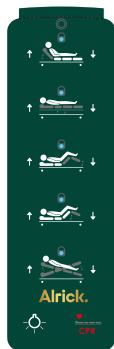

OPERATOR CONTROLS:

Electric Beds are controlled by a handset, which is plugged into the connection port. (Normally located below and at the edge of the mattress base towards the head end).

Main Bed Handset Code: HB8X016-00 (Second Bed Handset Code: J03508)



UPPER POSTURE, ELECTRIC BACKREST:

The top pair of push buttons control the electric upper posture as indicated on the handset.
Beware of potential pinch points when lowering backrest.



HI-LOW OPERATION:

The second pair of buttons on the handset control the smooth operation of raising & lowering the bed. When raising and lowering the bed, ensure that no person or part of his or her body is under the bed. Failure to do so could cause serious physical injury.



ELECTRIC KNEEBREAK:

The third pair of push buttons on the handset control the operation of the kneebreak. The calf section of the kneebreak is fitted with a prop. Lift the calf section and fold down the prop, engaging it in the channel provided at the foot end of the bed.
Beware of potential pinch points when lowering Kneebreak.



CHAIR POSITION:

The fourth set of buttons operate the chair function. This raises the backrest and kneebreak, then goes into reverse Trendelenburg (feet down).
If fitted with optional Electric Vascular, these buttons operate the electric vascular function and the Auto Chair Position is not available.



ELECTRIC TRENDLENBURG:

The last pair of push buttons on the handset control the operation of the tilt. Beware of potential pinch points when using the tilt function.



UNDER BED LIGHT:

The last button on the bottom left controls the under bed light.

CASTORS



STANDARD CASTORS:

Brakes should only be released when the bed is moved.
Red pedal castors: total brake.
Green pedal castors: directional.



CENTRAL LOCKING CASTORS:

Brakes should only be released when the bed is moved.
To brake: press lever down.
To direction: lift lever all the way up.
To free wheel: centralise lever.

PULL OUT EXTENSION



STEP 1
Pull down red knobs.



STEP 2
Pull out foot board.



STEP 3
Lock in red knobs.



STEP 4
Insert metal infill.

OPTIONS & EXTRAS



HEAD AND FOOTBOARDS:

The standard head & footboard frames push in from end, allowing for rapid removal in an emergency.
WARNING: Please ensure the power cord is not in the way when fitting or removing end boards.



CPR:

The CPR is only for emergency use.
Pull red lever to allow backrest to go down manually.



BATTERY BACK UP:

This allows you to safely change the patients position during a power shortage.

SAFETY & WARNINGS

- Electric Beds should ONLY be used where earth leakage is fitted to the room or building.
- If any cords are cut or damaged the bed MUST not be used until checked, fixed and tagged.
- No patient, service personnel or any person must lie underneath an Alrick bed while power is turned on, if not obeyed this may result in serious bodily injury or death if the bed controls are actuated.
- Please check no person or object is under the bed before using.
- Hand control to only be given to fully competent patients.
- Never leave hand control unattended with bed switched on unless patient is fully competent.
- Fit optional nurse lockout for safety of incompetent patients.
- Only one person is to be on the bed at any time.

MAINTENANCE

- Beds should be inspected monthly for abnormal wear. For example, bolts loose or missing, castors, rubbers or brakes damaged.
- Rails must be inspected regularly for any wear or movement that may cause injury, do not use if unsure.
- If one or more of the bed functions is found not working correctly, notify the manufacturer for advice.

EVERY 6 MONTHS:

1. Check there is no visible damage to electrical system or mechanical parts.
2. Check wiring is secure and not interfering with bed movement.
3. If actuator shafts are not clean - wipe chromed shaft of actuators to remove lint, grease & dirt.
4. Extended actuators fully. (Do not wash as this removes fine layer of oil).
5. Correct operation of all bed functions including castor locking.
6. Check all accessories including rails are still functioning properly.
7. Check all pivot points on the bed for wear. Replace pivot pins or bolts and removable bushes (where fitted) if a pivot has a significant play.
8. Check all functions still run smoothly.


ANNUALLY:

- Professional service of product to be carried out annually & documented in 'Record of Performed Service Procedures and Regular Service Maintenance in the back of this instruction manual.


CLEANING

- Cleaning can be done often - using domestic detergents.
- Do not use acids, alkaline or any other chemical solvents
- When cleaning actuators, disconnect from mains power before wiping down with a moist cloth.
- When cleaning, be careful of entrapment, as it could cause physical injury.

2300 Tilt partner bed series (BED 1Page 1)

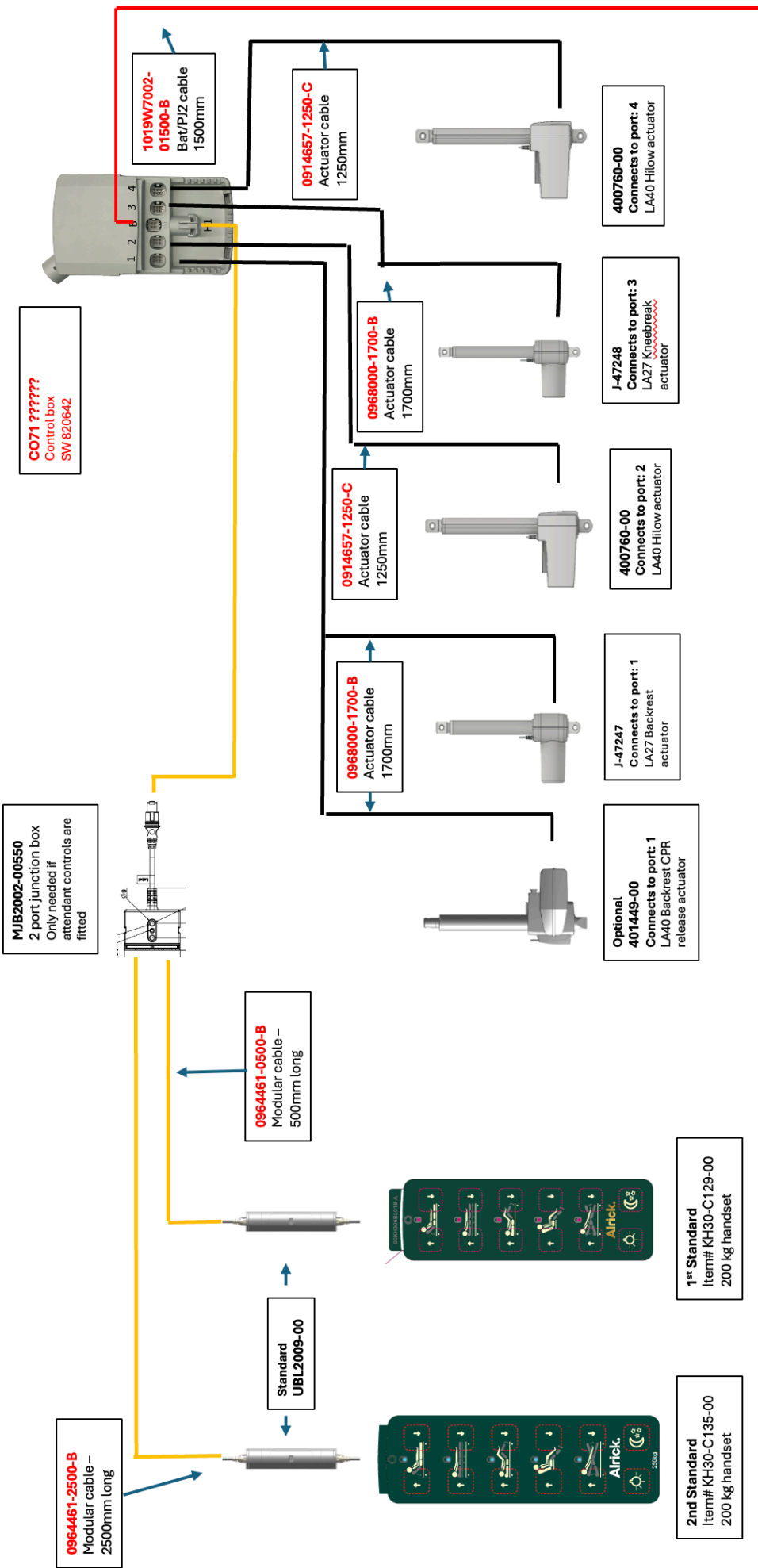

ELECTRICS (Bluetooth Handset - optional)

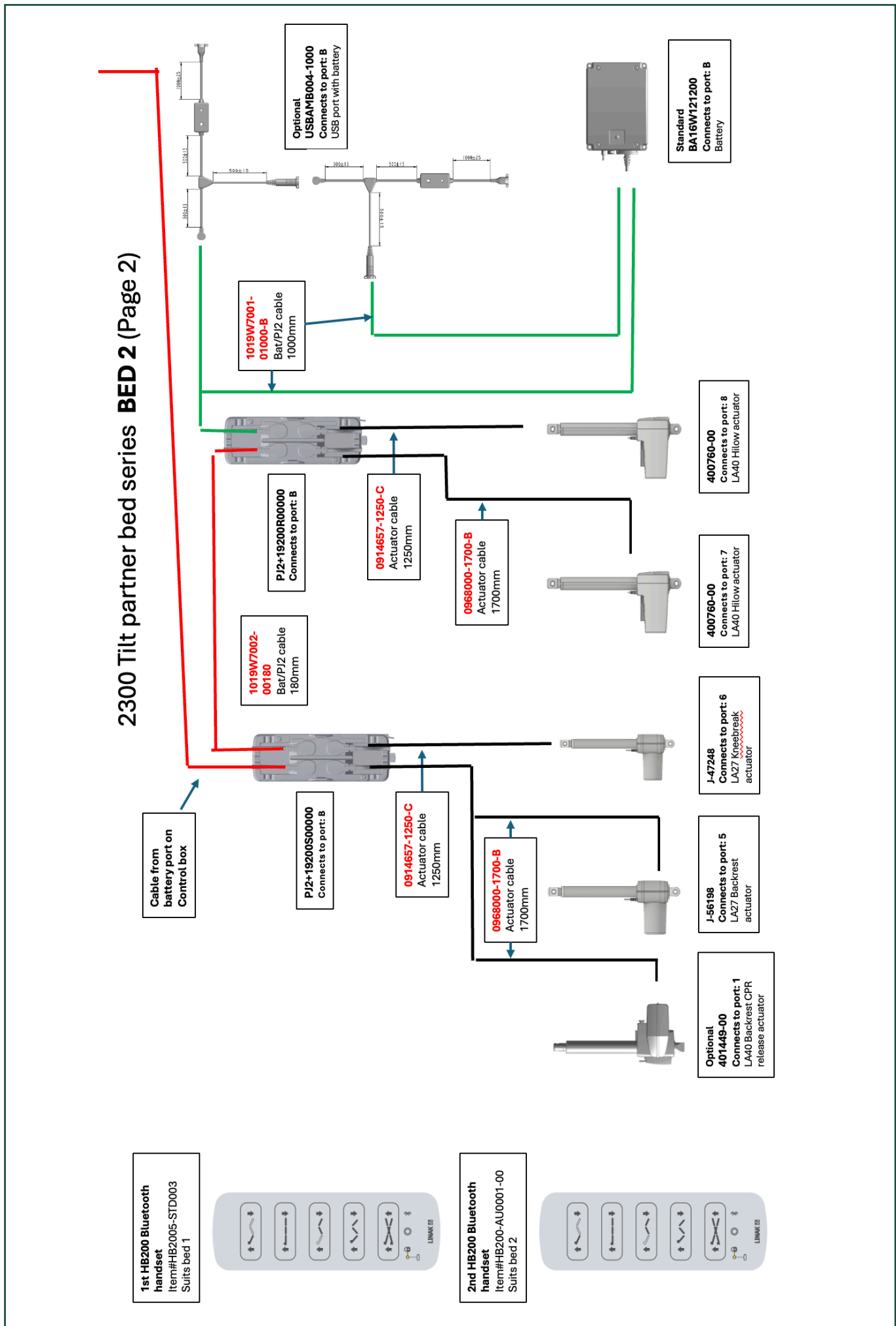

PARTNER BED ASSEMBLY

CONTENTS LIST

- 1 x Headboard
- 1 x Footboard
- 8 x Brackets (attached to bed)
- 16 x Sleeve nut bolt
- 16 x M6x16 button head screws
- 1 x Junction box cover
- 1 x Junction box clip
- 1 x Junction box plug


TOOLS

- 4 mm Allen Key




STEP 1

Remove the red 4 brackets from the head end of the bed.

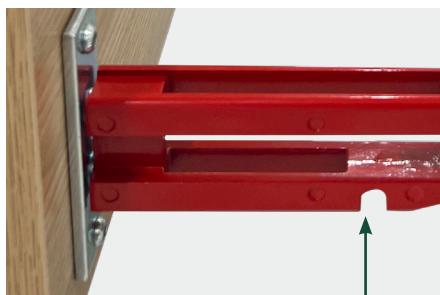

STEP 2

Lay headboard flat on the ground & attach the 4 brackets, removed in step 1 using the sleeve nut bolts and allen key. Ensure the attachment slot is facing down towards the bottom of the board.

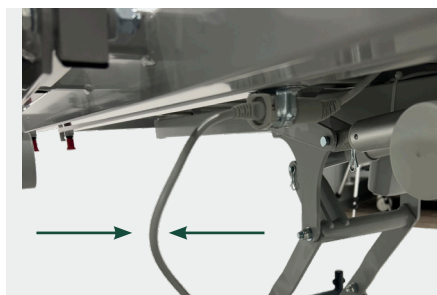

STEP 3

Position beds side by side in the highest position, with the 2 control cables at the centre of the two beds.



STEP 4

Insert the control cable into the junction box on the head frame of the bed.



STEP 5

Insert junction box plug (dust cap).



STEP 6

Clip junction box clip on.



STEP 7

Insert the functionality cable into the junction box under the back rest.



STEP 8

Clip on junction box cover.



STEP 9

Unlock locking knobs by twisting the 2 locking knobs at head end of each bed.



STEP 10

Slide red brackets attached to head board into the tube at the head end of the bed.



STEP 11

Twist 2 locking knobs on each bed, at head end to lock the head board into place.



STEP 12

Unlock locking knobs by twisting the 2 locking knobs at foot end of the bed.



STEP 13

Pull out brackets.



STEP 14

Hold end board in place and insert the top button head screw & sleeve nut, then tighten enough to hold in place.



STEP 15

Repeat step 14 on the other side of the end board.



STEP 16

Now the end board is in position, insert all screws and sleeve nuts. Don't tighten these fully until all nuts and bolts are in position.



STEP 17

Twist two locking knobs at foot end to lock the head board into place.



STEP 18

Insert extension infills if required.



Pictured without the
pull out extension.



Pictured with the
pull out extension.

LOCKING AND UNLOCKING LOCK OUT FUNCTIONS (optional)


HB400 HANDSET:

1. Locate the fob over the top of the handset as shown.
2. Then press the down key of the function to be locked out or unlocked. When locked, the LED function will be illuminated.




ACOM ATTENDANT CONTROL:

1. Press and hold the key button.

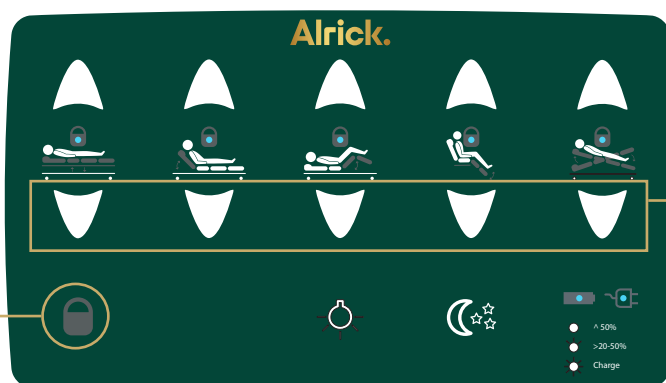

2. Then press the down key of the function to be locked or unlocked.

RESETTING THE BED FITTED WITH A CO71 LINAK CONTROL BOX:

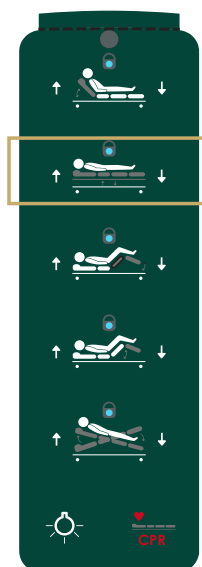


Press both hi-low buttons (up and down) together

- **It will beep for 5 to 10 seconds**
- don't let go of the buttons once beeping has stopped.
- **After another 10 seconds the bed will start to rise, it will go to the maximum height and return to the minimum height.**
- through this process don't let go of the buttons
- **Once at the lowest position the back and knee function motor will retract right in, once this is complete, you will get a double beep.**
- only let go of the buttons once you have heard these beeps


TROUBLESHOOTING

TROUBLESHOOTING WITH A HB400 HANDSET:

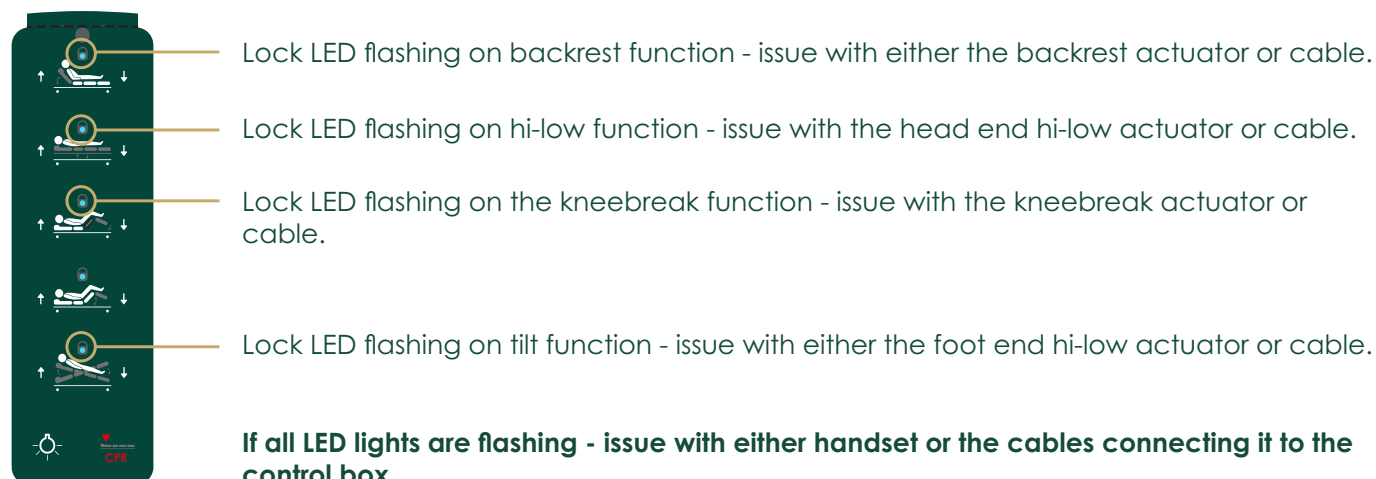

TROUBLESHOOTING WITH ACOM ATTENDANT CONTROL:

If all LED lights are flashing - issue with either handset or the cables connecting it to the control box.

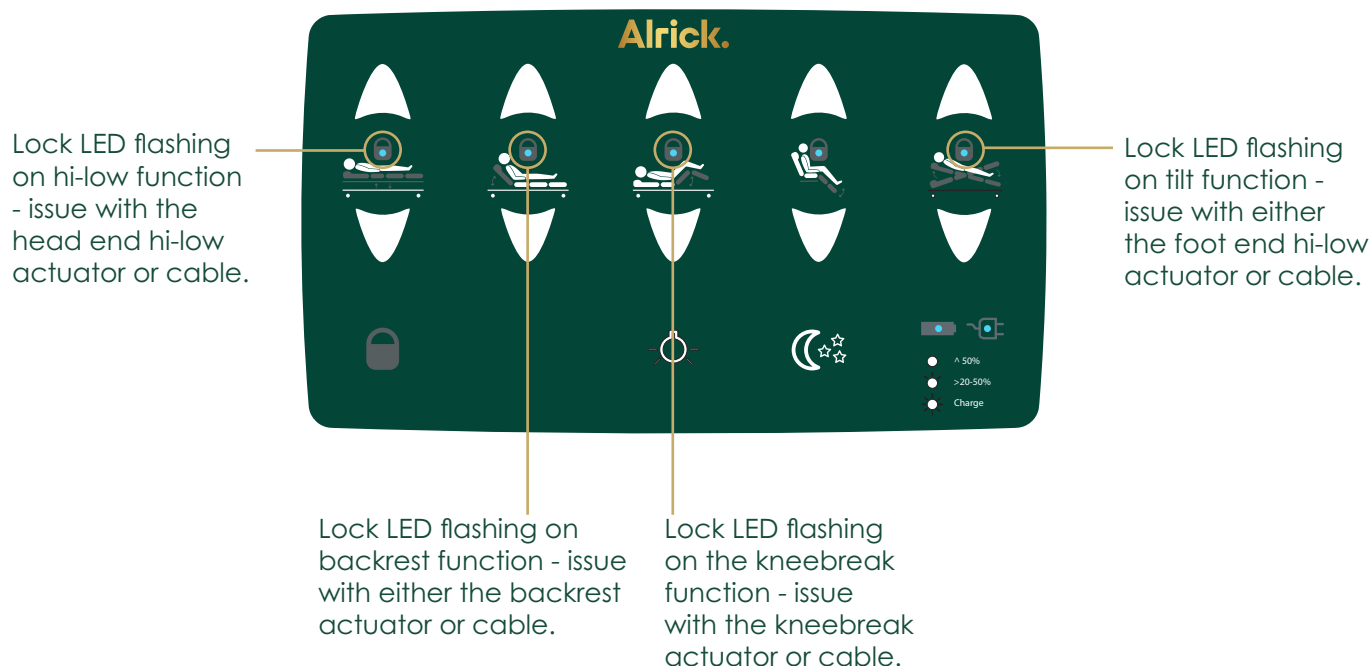


PAIRING BLUETOOTH DEVICES (optional)

DIRECT PAIRING

Direct Pairing is used for pairing a LINAK control directly to a LINAK control box that supports BLE.

1. Enter pairing mode
2. Move the hand control closer to the control box you want to pair it with
3. Pair the hand control with the control box

ENTER PAIRING MODE

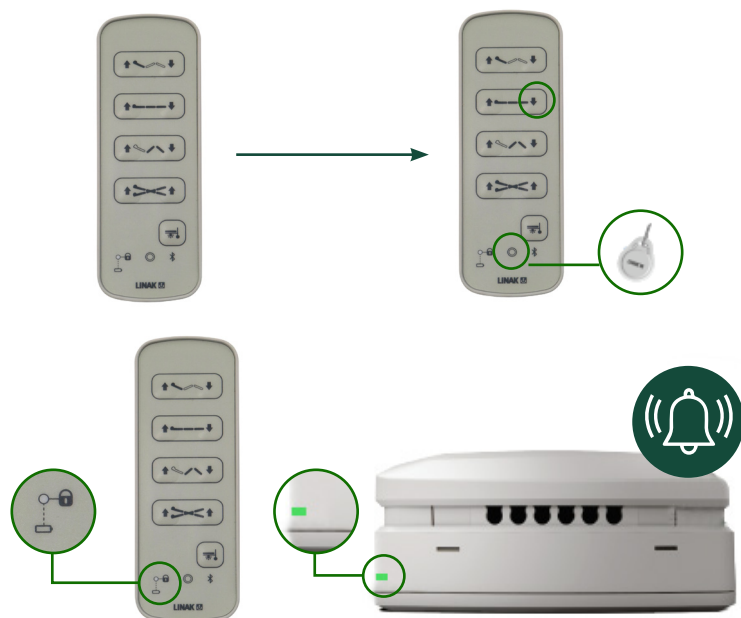


Place the magnet key over the **Hall Element (O)** and press **Hi-Lo Down Button** for 3 seconds to enter the pairing mode.

When in pairing mode, the control box buzzer will begin to beep, and the LED will begin to blink. The control unit LED will also begin to blink with the same frequency as the nearest control box.



In pairing mode, the light/sound frequency will increase when the HB200 gets closer to a control box.

CONNECTING TO THE CONTROL BOX

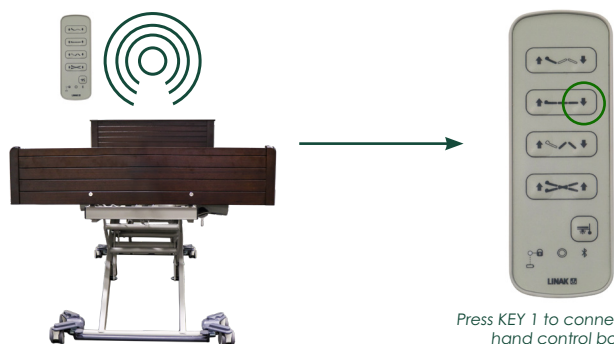

Press KEY 1 to connect the hand control box

When the hand control LED is blinking fast and the control box gives a high frequency sound in the same speed, the devices are ready for pairing. To finalise the pairing, press **Hi-Lo Down Button** on the hand control. Then, adjust your application with the handset.

RECORD OF PERFORMED SERVICE PROCEDURES AND REGULAR SERVICE MAINTENANCE

Date	Purpose & Description	Performance
------	-----------------------	-------------

Date	Purpose & Description	Performance
------	-----------------------	-------------

Date	Purpose & Description	Performance
------	-----------------------	-------------

RECORD OF PERFORMED SERVICE PROCEDURES AND REGULAR SERVICE MAINTENANCE

Date	Purpose & Description	Performance

I confirm that technical staff have undertaken the repairs as described in the Record of Performed Service Procedures and Repairs Maintenance, and the equipment is now functioning properly.

Date _____ Customer Signature _____

ALRICK 2300 PARTNER BED WARRANTY

CONDITIONS

Your Alrick 2300 Partner Bed is guaranteed for 10 years on the frame from the date of purchase, 5 years on the electrics (LINAK), 1 years on castors, and 1 year on the battery. Professional service of product to be carried out annually and documented in the 'Record of Performed Service Procedures and Regular Service Maintenance'. This warranty applies only on the parts manufactured or approved by Active Medical.

EXCLUSIONS

This warranty does not extend those items that may need replacement due to normal wear and tear, or to damage caused by misuse or accident, for which Active Medical or its distributors cannot be held responsible.

If your bed should need attention due to failure as defined under the warranty above, please contact the distributor from whom you purchased your bed.

Please fill in the below and retain for your own records.

Date Purchased

Model

Serial Number(s)

Distributor

Distributor Address

Distributor Contact Number

ACCESSORIES



Horizontal
400HR Deluxe Bed Rails
Warranty: 12 months



Clamp on Support Rail
Warranty: 12 months



Combo Grab Handle
Warranty: 12 months



Emergency Bed Power Unit
Warranty: 12 months



StandEzy Straight Handle
Warranty: 12 months



StandEzy Hook Handle
Warranty: 12 months



IV Pole
Warranty: 12 months



Foam Bolster
Warranty: 12 months



Self Help Pole
Warranty: 12 months



Bariatric Self Help Pole
Warranty: 12 months



T-Bar Buffer
Warranty: 12 months



Reading Lamp
Warranty: 12 months

Issue Date: 241024
Revision Date: 130826
V6

Alrick.

1800 282 001
hello@alrick.au
www.alrick.au