


Raz AT Tilt Adjustable Height Shower Commode

Product Code: 1948

Available with so many adjustments & options, the Raz-AT one of the easiest tilting mobile shower chairs to operate.

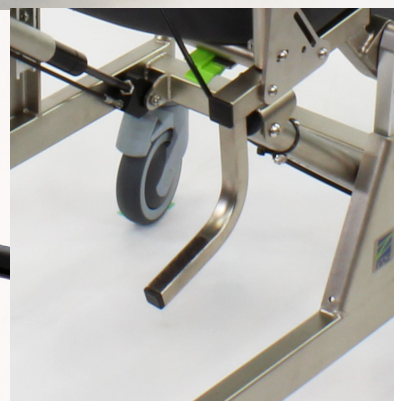

The Raz-AT (Attendant Tilt) mobile shower commode chair tilts up to 40° to provide extra pressure reduction and comfort. The unique frame design utilizes locking gas springs that help lift the seat to an upright position.

The Raz-AT mobile shower commode chair comes standard with the Ischial & Pelvic Alignment System (IPAS), 4" of height adjustment, adjustable-tension back upholstery, molded headrest and fixture, 5" dual-locking castors, flip-up padded armrests, and a commode pan. The easy-to-clean, polished stainless steel frame fits over most toilets and is designed for easy side and front access by the user or attendant.

Product Features:

- Adjustable seat height
- Flip Up Padded arm rests
- Adjustable angle backrest with adjustable tension
- 4 x Dual-locking (swivel and roll) castors
- Ischial & Pelvic Alignment System (IPAS), which allows the commode seat aperture to be comfortably centered around the ischial tuberosities.
- Adjustable / Quick-Release Seat
- Locking gas struts to facilitate tilting

Specifications:

Product Code	1948
Frame Width	460 mm
Overall Width	570 mm
Overall Length	990 mm
Seat Width	460 mm
Seat Depth	400 mm / 480 mm
Clearance Height (Under Seat) (with / without pan brackets)	430 - 530 mm 450 - 640 mm
Armrest Height (From Seat)	190 - 270 mm
Seat to Floor Height	540 - 640 mm
Product Weight	28 kg
SWL	160 kg