

Secura 3 Mental Health Bed

Product Code: 7-SEC3

The ultimate safe and secure solution, eliminating access to the patented bed mechanisms, with the bed surround maintaining the concealment in both, Hi and Low positions, the lower bed section sitting stationary with the bed. The sides of the bed are recessed to allow for nurses to kneel on the bed without the obstruction in the shin area.



Product Features:

- Australian Owned & Made
- Electric Hi-Low
- Enclosed anti-ligature area
- Embedded footboard control
- Concealed wiring
- Enclosed mattress supports
- Global Lock Out with Fob
- Braided electrical power supply
- 75mm Central Locking Castors
- Concealed Castors
- Headboard & footboard lower than mattress height



Brake Lever
Access



Anti Tamper
Screws



13mm Natural Oak
Compact Laminate
Surround



Ironstone
Powdercoat

Specifications:

Size	Single	King Single
Overall Dimensions	930 x 2030 mm	1100 x 2030 mm
Sleeping Surface	900 x 2000 mm	1070 x 2000 mm
Low Position	400 mm	400 mm
High Position	700 mm	700 mm
SWL	250 kg	250 kg

Issue Date: 010524
Revision Date: 290525 | V4



Backrest is optional

Ensuring high comfort & safety in mental health environments.

Optional Features - Added at time of manufacture

- Kneebreak
- Backrest
- Braided electrical power supply with fiberglass shroud
- Built in 200mm extension
- Battery back up
- Custom powdercoat colours
- Custom surround colours
- Extra length
- Removable Mains Cable option (under Key Lock Out)


Optional Mobile Battery Trolley

- Added at time of manufacture


- Providing a means to lower and raise a bed when a room is required to have no power.
- For use in high risk facilities where rooms don't have power outlets
- On wheels making it easily mobile
- Able to be mounted onto a linen trolley
- Raises bed to allow access to brake lever
- Raises bed to prevent staff injury when caring for a patient or making the bed

Surround Colours

Custom Compact Laminate Surround Colours Available

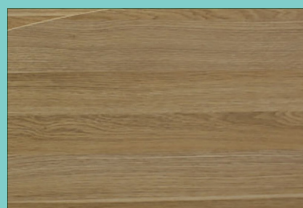

Natural Oak



Oyster Grey

Accessories - Added after manufacture

- Mattresses
- Bluetooth controls

