

Staebel Forte 320 Patient Lifter

Product Code: 2898

The ultimate heavy duty patient lifting hoist for all healthcare environments.



The lifter of choice for institutional use where a high capacity bariatric lifter is required.

The Staebel Forte 320 Patient Lifter is a compact, ergonomic designed heavy-duty mobile patient lifting hoist with a safe working load of 320kg. With robust powder coated steel construction this lifter is ideally suited for institutional use or as a compact heavy-duty home care hoist for the bariatric patient. The Forte 320, whilst being compact, offers superior patient lift space through the combination of the recessed actuator and the use of SPACE LINK™ technology. The unique feature of SPACE LINK™ allows a patient to be rotated through 360° at full lift height, thus making the carers job easier and safer.

The Forte 320 Lifter provides infection control advantages through the aesthetic design of the concealed electric leg spread mechanism. The Forte 320 lifting range allows patients to be lifted from the floor, bed or chair with ease and safety. The Forte 320 is a powered unit using trusted industry leading LINAK components, giving a safe working load (SWL) of 320kg.

The Staebel Forte 320 Patient Lifter is the lifter of choice where a high capacity bariatric lifter is required.

Product Features:

- Supplied complete with battery, in-built battery charger and instruction manual.
- Safe Working load of 320kg
- Battery life indicator on lifter as well as handset
- Additional up and down buttons on control box
- Australian Owned and Made
- Compatible with our Staebel Slings & a range of hoist attachments



Standard Yoke



4 Point Yoke




Tri Hook Yoke



Pivot Frame



Powered Pivot Frame

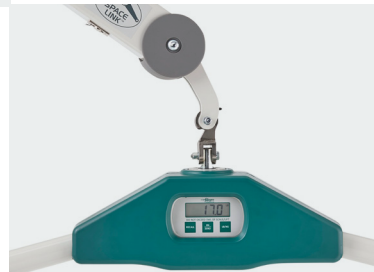

Integrated Weigh Scale

Specifications:

Product Code	2898	Under Bed Clearance	111 mm
Base Length	1420 mm	Minimum Yoke Height	245 mm
External Base Width Closed	715 mm	Maximum Yoke Height	1940 mm
External Base Width Open	1200 mm	Lifter Weight	63 kg
Internal Base Width Closed	600 mm	Minimum Storage Height	1445 mm
Internal Base Width Open	1080 mm	SWL	320 kg