

Relax Easy Cushion

Product Code: 351REL

The ideal pressure relieving cushion for use in the healthcare environments.



Ideal pressure relieving cushion for use in nursing home environments.

As the name suggests, this cushion from the Relax range is easy to use and is an ideal pressure relieving cushion for use in the nursing home environment. It has no designated up or down, or front or back, so can simply be placed into position. It features a layer of FloGel sandwiched between an upper and lower layer of viscoelastic foam. This combination of pressure relieving FloGel and viscoelastic foam facilitates effective pressure re-distribution while allowing the patient to immerse themselves comfortably into the cushion. It conforms to the shape of the user, increasing surface area and relieving pressure.

Product Features:

- 1.9" (50mm) thick
- Suitable for medium – high risk clients
- Clear protective PU film coating is easy to wipe clean and is suitable for incontinent users
- Lightweight and easy to transport
- Fully rotational design - it has no top, bottom, front or back
- Includes a black waterproof wipeable PU stretch outer cover

P



Prevention
Cushion




Combination of
FloGel & Foam



Medium – High Risk
Category

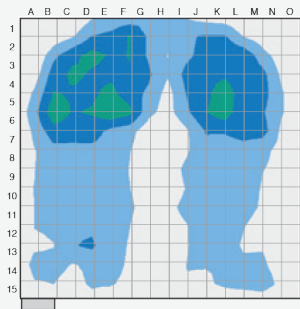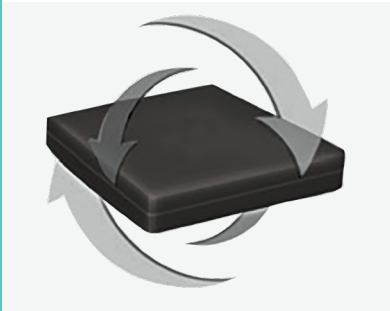

Pressure diagram with Relax
Easy cushion

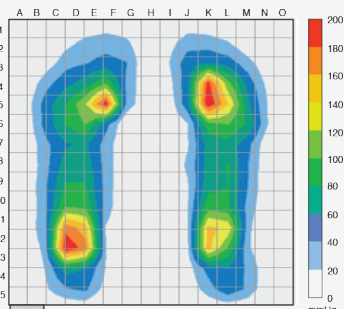

Pressure diagram without
Relax Easy cushion

Specifications:

Product Code	35100	35312	35313	35315	353120
Width	400 mm	450 mm	450 mm	500 mm	500 mm
Length	400 mm	400 mm	450 mm	450 mm	500 mm
Thickness	50 mm	50 mm	50 mm	50 mm	50 mm
SWL	130 kg	130 kg	130 kg	130 kg	130 kg