

# Staebel Relax DuoGel Cushion

Product Code: 353RDG

The ideal pressure relieving cushion for long-term use in all healthcare environments.



# The Relax Duogel effectively provides pressure relief & weight distribution.

The Relax Duogel is a unique combination of TruGel and FloGel working together effectively to provide excellent pressure relief and weight distribution. The Duogel’s anatomical pre-contoured base provides for medial and lateral thigh positioning as well as support for the greater trochanter of the thigh bone. The thigh troughs are lined with TruGel and the ischial area is lined with FloGel which is set on a viscoelastic foam giving greater comfort and seating stability.

## Product Features:

- 2.9” (75mm) thick
- The anatomical pre-contoured base provides for medial and lateral thigh positioning as well as support for the greater trochanter of the thigh bone
- The thigh troughs are lined with TruGel and the ischial area is lined with FloGel which is set on a viscoelastic foam giving greater comfort and seating stability
- Clear protective (PU) film coating with no seams is easy to wipe clean and is suitable for incontinent users
- Suitable for those high risk of pressure injury.
- Includes a black waterproof wipe down outer cover
- Safe Working Load of 130kg

P


Prevention Cushion




Combination of TruGel & FloGel



High Risk Category

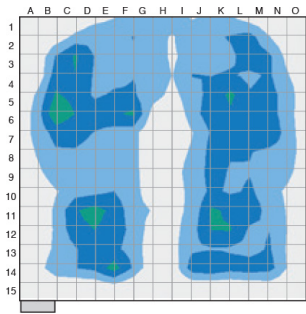

Pressure diagram with Relax Duogel cushion

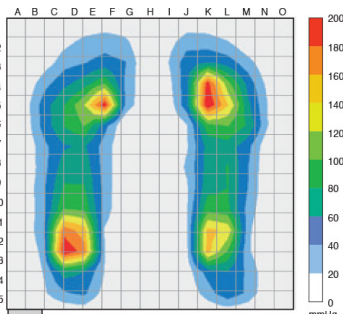

Pressure diagram without Relax Duogel cushion

## Specifications:

Product Code	35300	35302	35303	35304	35307
Width	400 mm	450 mm	450 mm	500 mm	500 mm
Length	400 mm	400 mm	450 mm	450 mm	500 mm
Thickness	80 mm	80 mm	80 mm	80 mm	80 mm
SWL	130 kg	130 kg	130 kg	130 kg	130 kg