

Staebel Relax DuoGel Cushion

Product Code: 353RDG

The ideal pressure relieving cushion for long-term use in all healthcare environments.



130kg
SWL



High Risk
Cushion



Prevention
Cushion



Combination Gels
for Target Pressure Relief




Pre-Contoured
Surface

The Staebel Relax DuoGel Cushion features a unique combination of TruGel and FloGel technologies, designed to deliver both effective pressure redistribution and enhanced seating stability. FloGel, positioned in the posterior region, facilitates immersion and envelopment around the ischial tuberosities and sacral area, supporting pressure relief over bony prominences. TruGel, integrated into the medial and lateral thigh troughs, promotes pelvic stability and controlled positioning of the lower limbs. This dual-gel configuration, set atop a viscoelastic foam base, offers targeted support for the greater trochanters and contributes to overall seating comfort, postural alignment, and skin integrity preservation.

To preserve foam integrity, the inner polyurethane (PU) film coating provides an additional barrier against fluid ingress, enhancing the longevity of the memory foam, particularly beneficial for users with incontinence.

Product Features:

- 2.9" (75mm) thick
- Dual-layer gel construction featuring FloGel and TruGel for combined pressure relief and stability
- FloGel in the posterior region allows immersion and envelopment over ischial tuberosities and sacrum
- TruGel in the thigh troughs enhances pelvic stability and controlled leg positioning
- Supports medial and lateral thigh alignment and offloads the greater trochanters
- Viscoelastic foam base for improved comfort and sustained pressure redistribution
- Designed to promote postural alignment and help maintain skin integrity
- Clear protective (PU) film coating for incontinent users
- Suitable for those high risk of pressure injury
- Outer cover can be machine washed (60°C)
- Includes a black, wipeable, 4-way stretch, vapour permeable, water resistant outer cover
- WARRANTY: 2 Years

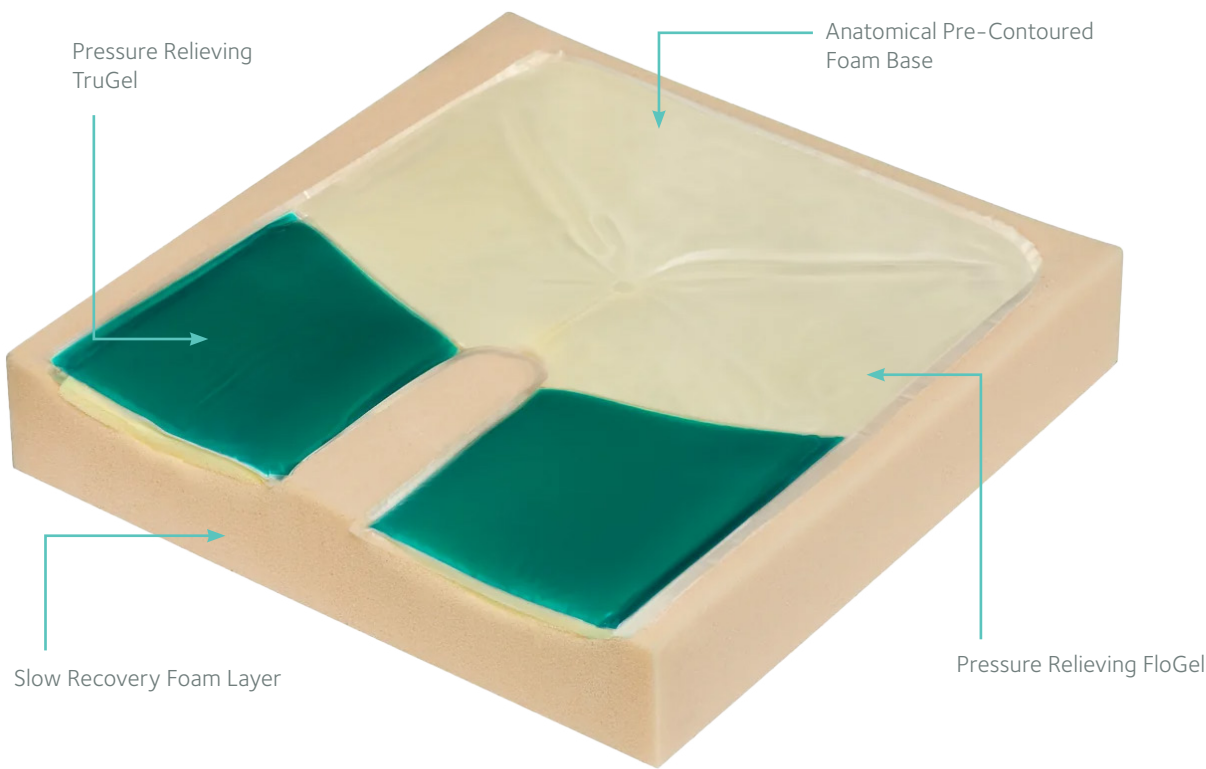

Specifications:

Product Code	35300	35302	35303	35304	35307
Width	400 mm (16")	450 mm (18")	450 mm (18")	500 mm (20")	500 mm (20")
Length	400 mm (16")	400 mm (16")	450 mm (18")	450 mm (18")	500 mm (20")
Thickness	75 mm (2.9")	75 mm (2.9")	75 mm (2.9")	75 mm (2.9")	75 mm (2.9")
SWL	130 kg	130 kg	130 kg	130 kg	130 kg