

Staebel Relax Easy High Profile Cushion

Product Code: 353REH

The ideal pressure relieving cushion for long-term use in all healthcare environments.






The ideal pressure relieving cushion for long term use in all healthcare environments.


This cushion from the Relax range is easy to use and is an ideal pressure relieving cushion for use in a long term care environment. It features a layer of FloGel sandwiched between an upper and lower layer of viscoelastic foam. Relax High Profile has the same properties as the standard lax, therefore it conforms to the shape of the user, increasing surface area and relieving pressure.

Product Features:

- Fully Rotational – No top, bottom, front or back
- High density viscoelastic foam with silicone fluid gel at its core, provides superior pressure redistribution and immersion
- Silicone oil is centered between the two layers of slow foam, it contains preformed channels to make sure the oil is evenly distributed
- Thanks to additional height and materials, the Easy High profile offers pressure relief for up to 190kg
- Conforms to the shape of the user, increasing surface area and relieving pressure
- Increased depth of memory foam and silicone oil allows for greater immersion and comfort
- Combination of materials reduces 'hammock' effect
- Ideal for a patient requiring additional stability, and/or greater cushion depth in their wheelchair

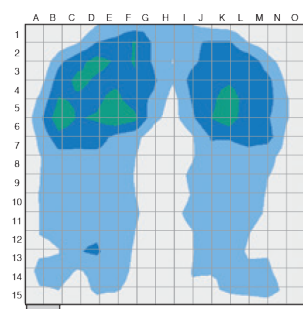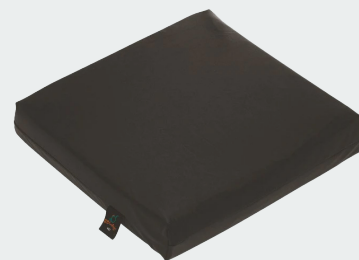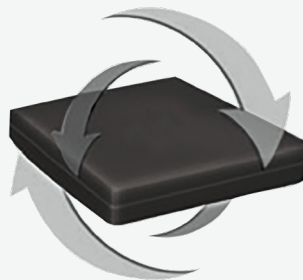

Pressure diagram with Relax Easy cushion

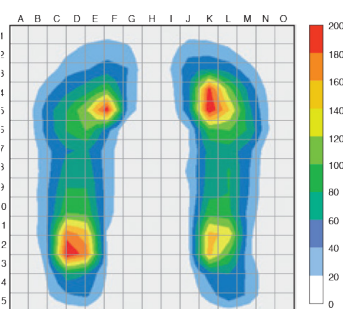

Pressure diagram without Relax Easy cushion



Prevention
Cushion



Combination of
FloGel & Foam



Medium Risk
Category

Specifications:

Product Code	35392	35394	35395	35396	35397
Width	400 mm	450 mm	450 mm	500 mm	500 mm
Length	400 mm	400 mm	450 mm	450 mm	500 mm
Thickness	85 mm	85 mm	85 mm	85 mm	85 mm
SWL	190 kg	190 kg	190 kg	190 kg	190 kg