

Staebel Relax Easy High Profile Cushion

Product Code: 353REH

Staebel Relax Easy High Profile Cushion designed for the aged and home care settings, featuring a fully rotatable design, mitigating the risk of improper application and supporting consistent clinical outcomes.



190kg
SWL



At Risk - High Risk
Cushion



Prevention
Cushion



Fully Rotatable
Design



Aged & Home
Care Appropriate


The Staebel Relax Easy High Profile Cushion has a layer of FloGel between an upper and lower layer of viscoelastic foam. The combination of pressure-relieving FloGel and viscoelastic foam facilitates effective pressure redistribution by immersion and envelopment, conforming to the shape of the user and increasing surface area.

The Relax Easy High Profile is well suited for use in the aged care and home care settings as it features a fully rotatable design, mitigating the risk of improper application and supporting consistent clinical outcomes. To preserve foam integrity, the inner polyurethane (PU) film coating provides an additional barrier against fluid ingress, enhancing the longevity of the memory foam, particularly beneficial for users with incontinence.


Product Features:

- 3.3" (85mm) thick
- Suitable for at risk - high risk users
- Clear protective PU film coating is easy to wipe clean and is suitable for incontinent users
- Lightweight and easy to transport
- Outer cover can be machine washed (60°C)
- Includes a black wipeable, 4-way stretch, vapour permeable, water resistant outer cover
- WARRANTY: 2 Years

**Relax Easy Low Profile option available*

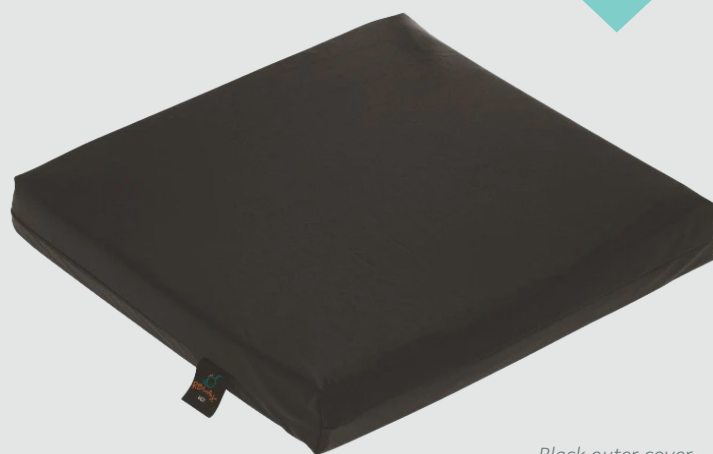

Black outer cover

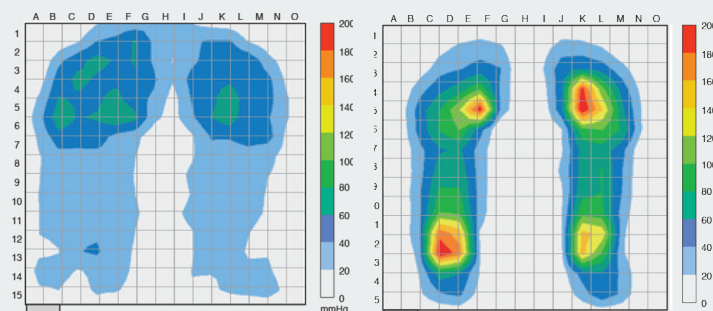

*Pressure diagram with
Relax Easy cushion*

*Pressure diagram without
Relax Easy cushion*

Specifications:

Product Code	35392	35394	35395	35396	35397
Width	400 mm (16")	450 mm (18")	450 mm (18")	500 mm (20")	500 mm (20")
Length	400 mm (16")	400 mm (16")	450 mm (18")	450 mm (18")	500 mm (20")
Thickness	85 mm (3.3")	85 mm (3.3")	85 mm (3.3")	85 mm (3.3")	85 mm (3.3")
SWL	190 kg	190 kg	190 kg	190 kg	190 kg