

Staebel Relax GelCell Cushion

Product Code: 353RGC

The Staebel Relax GelCell is a unique combination of AIR and GEL technology, enabling personalised pressure redistribution and comfort to the very high-risk user.



130kg
SWL



Very High Risk
Cushion



Prevention
Cushion



Microclimate
Management




Personalised
Immersion & Envelopment

The Staebel Relax GelCell Cushion is a unique combination of AIR and GEL technology, enabling personalised pressure redistribution and comfort to the very high-risk user. The Gelcell consists of 34 interconnected air-filled cells, each covered with a top layer of pressure-relieving silicone-based TruGel.

The immersion by use of air as the medium in the cushion maintains capillary function in the tissue, ensuring blood flow and oxygen uptake. When pressure is applied, air circulates to each of the cells, distributing weight evenly, thus reducing high-pressure concentrations around bony prominences. The GelCell has a positive effect on temperature management due to the interconnected air cell design, permitting air circulation across the cushion.

Product Features:

- 2" (50mm) thick
- Personalised immersion and envelopment
- Suitable for very high-risk users
- Lightweight and easy to transport
- Easy to inflate and adjust
- Outer cover can be machine washed (60°C)
- Includes a black, wipeable, 4-way stretch, vapour permeable, water resistant outer cover.
- Pump included
- WARRANTY: 2 Years

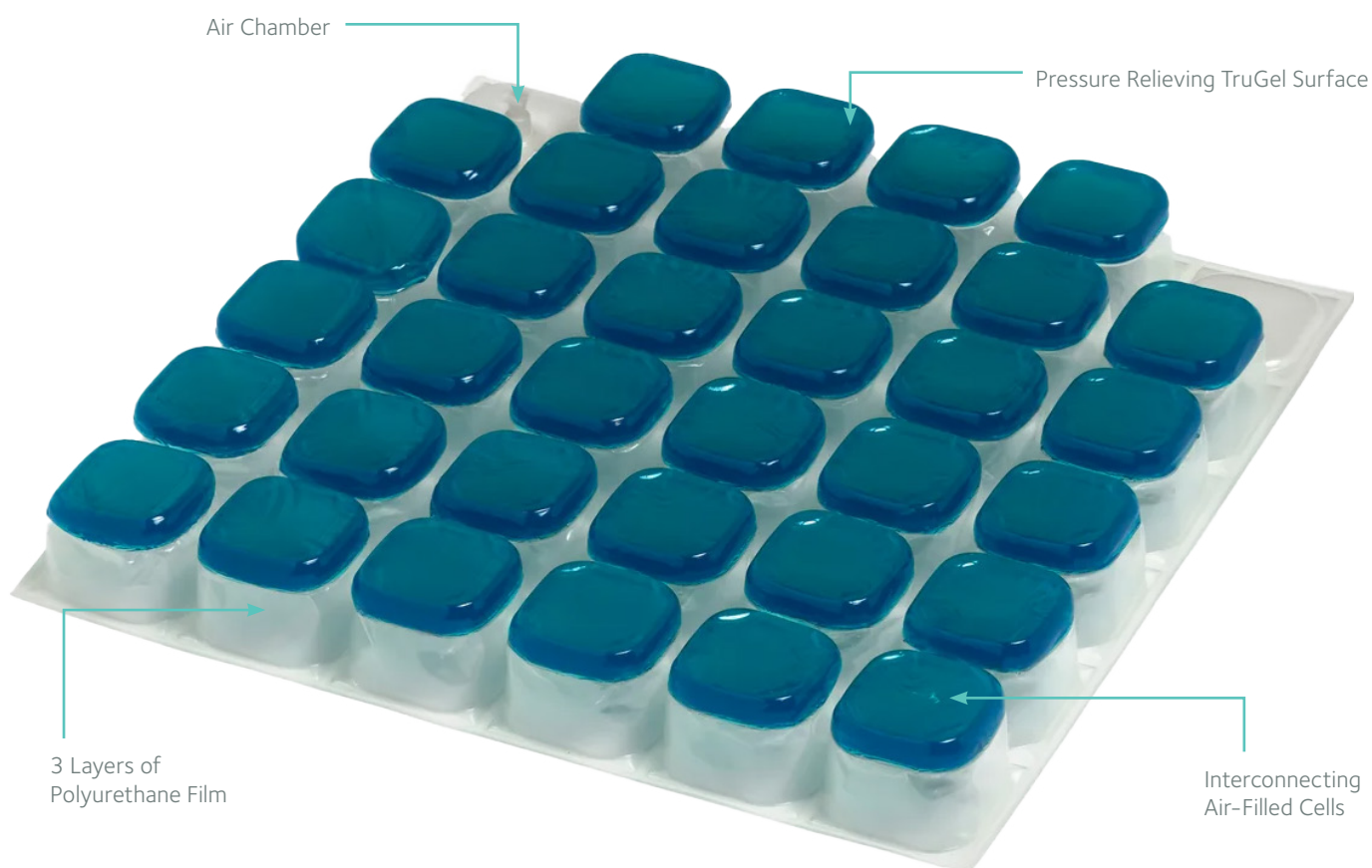

Specifications:

Product Code	35320	35321	35322	35323	3532320
Width	400 mm (16")	430 mm (17")	450 mm (18")	450 mm (18")	500 mm (20")
Length	400 mm (16")	430 mm (17")	400 mm (16")	450 mm (18")	500 mm (20")
Thickness	50 mm (2")	50 mm (2")	50 mm (2")	50 mm (2")	50 mm (2")
SWL	130 kg	130 kg	130 kg	130 kg	130 kg

Issue Date: 140425
Revision Date: 220725 | V2