

Staebel Relax Trio Cushion

Product Code: 353RTC

An ultra-lightweight wheelchair cushion that is suitable for those low to at-risk of pressure injuries.



80kg
SWL



Low - At Risk
Cushion



Prevention
Cushion



Ultra Lightweight
Design




Aged & Home
Care Appropriate

The Staebel Relax Trio is an ultralight cushion that benefits from four sections of FloGel being positioned into a polyurethane base, with a viscoelastic foam upper layer. It is specially designed to conform to the shape of the user, increasing surface area and relieving pressure. Suitable for low to at risk users, the low profile cushion can be rotated either up or down; it has no designated front or back.

To preserve foam integrity, the inner polyurethane (PU) film coating provides an additional barrier against fluid ingress, enhancing the longevity of the memory foam, particularly beneficial for users with incontinence.

Product Features:

- 1.5" (40mm) thick
- Recovery foam upper layer with a fast recovery foam core and a viscoelastic slow foam base.
- Rotatable cushion with no front or back.
- Clear protective (PU) film coating with no seams is easy to wipe clean and is suitable for incontinent users.
- Highly durable and lightweight.
- Suitable for those low to at-risk of pressure injury
- Outer cover can be machine washed (60°C)
- Includes a black wipeable, 4-way stretch, vapour permeable, water resistant outer cover
- WARRANTY: 2 Years

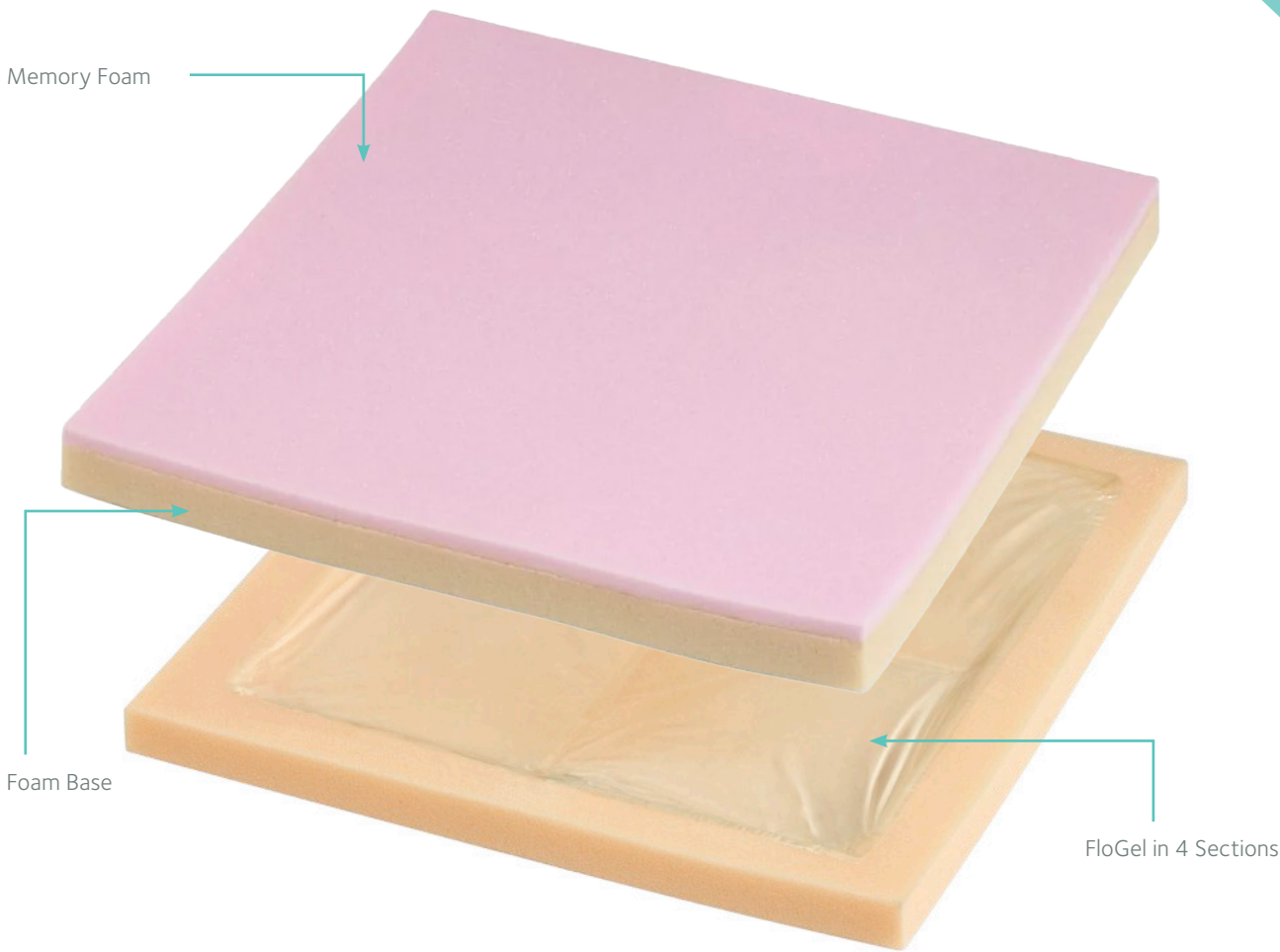

Specifications:

Product Code	35329	35330	35331	35332	35333
Width	360 mm (14")	400 mm (16")	430 mm (17")	450 mm (18")	450 mm (18")
Length	400 mm (16")	400 mm (16")	430 mm (17")	400 mm (16")	450 mm (18")
Thickness	40 mm (1.5")	40 mm (1.5")	40 mm (1.5")	40 mm (1.5")	40 mm (1.5")
SWL	80 kg	80 kg	80 kg	80 kg	80 kg

Issue Date: 140425
Revision Date: 220725 | V2