

Staebe! Total Flow Max Active Air Mattress

Product Code: 3-CM-M-FAS

A full replacement alternating mattress system, designed as a cost-efficient solution across various healthcare environments.



Product Features:

- Nylon TPU independently removable cells
- Cell-on-cell design to prevent bottoming out and increase stability when the alternating mode is active
- Remote control to enable independence
- Pillow function/static head cells
- CPR quick release tag & anchoring straps
- Multi-stretch, antimicrobial, machine washable, waterproof & vapour permeable PU cover
- Micro-air loss cells in the torso region to support microclimate management
- Transport Mode
- Warranty: Mattress - 12 months, Pump - 24 months
- Standard: ISO 9001:2015 | ISO 13485:2016

Pump Features:

- Cycle Time: 10, 15, 20, 25 minutes
- Control panel lock to prevent unauthorised or accidental tampering
- 20-minute automatic timer for max inflate mode to reduce risks of prolonged max inflation
- Memory function in the case of a power outage
- Air Output: 8 litres/minute
- Power Input: AC220-240V, 50-60Hz
- Power Consumption: 7W
- Fitted with hanging brackets
- Dimensions (WxLxH): 240 x 345 x 125 mm



30-200kg
SWL



High Risk
Mattress



Prevention
Mattress




Removable
Heel Cells



Micro Air
Loss

The Staebel Total Flow Max is a full replacement alternating mattress system, designed as a cost-efficient solution across various healthcare environments. The Staebel Total Flow Max has a cell-on-cell design with static base cells for stability and a static pillow region. The alternating air mattress is effective for the prevention of patients at high risk of developing pressure injuries. The cells in the heel region can be easily decoupled if personalised offloading is required for the individual. The design additionally supports the microclimate by having micro air-loss cells within the torso region.

With a simple-to-use digital design, this mattress allows the carer to quickly set up the mattress system and adjust it to the optimal pressure based on the user's specific weight. The Staebel Total Flow Max incorporates cutting-edge technology with constant auto-pressure sensing through an air sensor and microcomputer regulator to provide ultra-low interface pressure, promoting an even ambience for natural sleep patterns.

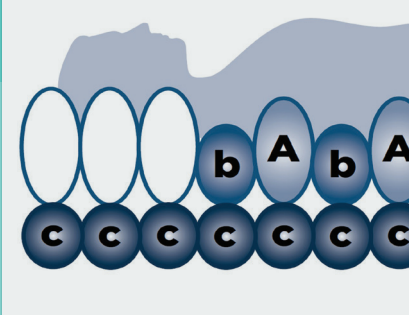

Specifications:

Product Code	3-CM-M-FAS-S	3C-M-M-FAS-KS
Size	Single	King Single
Width	900 mm	1070 mm
Length	2000 mm	2000 mm
Thickness	200 mm	200 mm
SWL	30 - 200 kg	30 - 200 kg

Control Panel:

