

Staebel OptiCell SMART R6 Active Air Mattress

Product Code: 327109SR6

The Staebel OptiCell SMART R6 Active Air Mattress sets out to improve the efficiency of the care staff and the time that has been spent on developing an easy-to-use and intuitive panel.



230kg
SWL



Very High Risk
Mattress

4

Suitable for
Stage 4

T

Treatment
Mattress



Micro Air
Loss



Personalised
Offloading

The Staebel OptiCell SMART R6 Mattress is designed with easy and quick individual cell decoupling, without affecting the mattress's function, as each individual cell has a built-in check valve. With this construct, the Staebel Opticell SMART R6 facilitates personalised and targeted offloading, whereby offloading is particularly relevant when actively treating pressure injuries. The Staebel Opticell SMART R6 mattress can aid in the prevention and treatment of pressure injuries up to and including stage 4.

The mattress offers a patented TrueAlternating mode, which means that the deflated cell rises to complete inflation before the next cell cycle occurs; this ensures no movement up or down. This is unique to the OptiCell range and reduces/ eliminates the end-user feeling nausea and motion sickness, a problem that can sometimes arise in traditional alternating mattresses. The fully automatic SMART pump improves efficiency and reduces human error with an easy-to-use and intuitive panel, eliminating the need for manual weight entry.




Product Features:

- The mattress system provides comfort, support and pressure redistribution through immersion and envelopment
- Personalised offloading can be achieved by cell decoupling anywhere across the support surface
- Micro Air Loss air cells provide good ventilation and help reduce heat and moisture, supporting the microclimate
- The OptiTex™ mattress cover is water resistant, vapour permeable, hydrostatic and fire rated
- Quality product from Europe
- The Staebel OptiCell SMART R6 Mattress is 6" high when max inflated; the OptiCell SMART R8 is 8" high

SMART Pump Features:

- Automatic adjustment to the user's weight and bed positioning
- Quiet and vibration-free
- Automatic 20-minute timer after max inflate/care mode is enabled
- Battery saves pump settings in case of power failure
- The user-friendly display has a layout that is easy to understand
- The pump has a smooth surface without cavities to facilitate daily cleaning and disinfection
- Performance data supporting immersion and envelopment capabilities
- WARRANTY: 2 Years

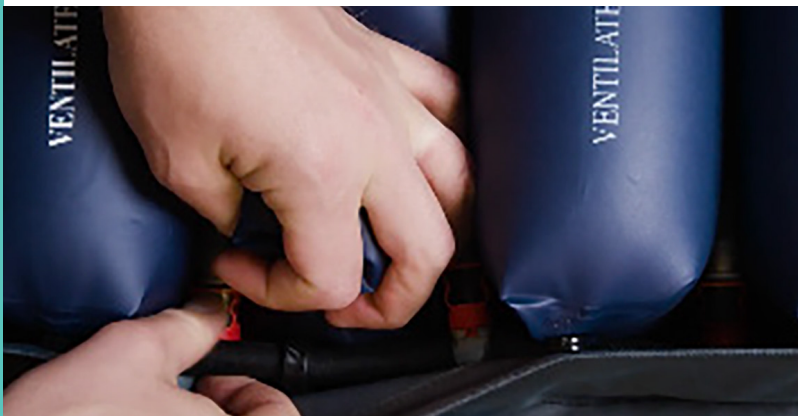

Specifications:

Product Code	3271091	3271092	3271093	3271094
Size	Single	Single	King Single	Euro
Width	850 mm	900 mm	1050 mm	1200 mm
Length	2000 mm	2000 mm	2000 mm	2000 mm
Thickness	160 mm	160 mm	160 mm	160 mm
SWL	230 kg	230 kg	230 kg	230 kg