

Staebel OptiCell SMART R8 Active Air Mattress

Produce Code: 327111SR8

The Staebel OptiCell SMART R8 Active Air Mattress sets out to improve the efficiency of the care staff and the time that has been spent on developing an easy-to-use and intuitive panel.



250kg
SWL



Very High Risk
Mattress

4

Aids in the Treatment of up
to & including **Stage 4 PI**



Micro Air
Loss



Personalised
Offloading

The Staebel OptiCell SMART R8 Mattress is designed with easy and quick individual cell decoupling, without affecting the mattress's function, as each individual cell has a built-in check valve. With this construct, the Staebel Opticell SMART R8 facilitates personalised and targeted offloading, whereby offloading is particularly relevant when actively treating pressure injuries. The Staebel Opticell SMART R8 mattress can aid in the prevention and treatment of pressure injuries up to and including stage 4.

The mattress offers a patented TrueAlternating mode, which means that the deflated cell rises to complete inflation before the next cell cycle occurs; this ensures no movement up or down. This is unique to the OptiCell range and reduces/eliminates the end-user feeling nausea and motion sickness, a problem that can sometimes arise in traditional alternating mattresses. The fully automatic SMART pump improves efficiency and reduces human error with an easy-to-use and intuitive panel, eliminating the need for manual weight entry.




Product Features:

- The mattress system provides comfort, support and pressure redistribution through immersion and envelopment
- Personalised offloading can be achieved by cell decoupling anywhere across the support surface
- Micro Air Loss air cells provide good ventilation and help reduce heat and moisture, supporting the microclimate
- The OptiTex™ mattress cover is water resistant, vapour permeable, hydorstatic and fire rated
- Quality product from Europe
- The Staebel OptiCell SMART R6 Mattress is 6" high when max inflated; the OptiCell SMART R8 is 8" high

SMART Pump Features:

- Automatic adjustment to the user's weight and bed positioning
- Quiet and vibration-free
- Automatic 20-minute timer after max inflate/care mode is enabled
- Battery saves pump settings in case of power failure
- The user-friendly display has a layout that is easy to understand
- The pump has a smooth surface without cavities to facilitate daily cleaning and disinfection
- Performance data supporting immersion and envelopment capabilities
- WARRANTY: 2 Years

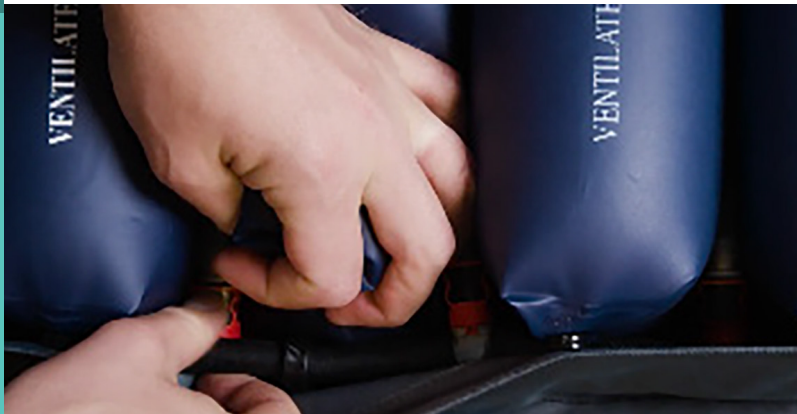

Specifications:

Product Code	Size	Width	Length	Thickness	SWL
3271111	Single	850 mm	2000 mm	200 mm	250 kg
3231112	Single	900 mm	2000 mm	200 mm	250 kg
3231113	King Single	1050 mm	2000 mm	200 mm	250 kg
3271114	Euro	1200 mm	2000 mm	200 mm	250 kg
3271118	Long Single	850 mm	2100 mm	200 mm	250 kg
3271119	Long Single	900 mm	2100 mm	200 mm	250 kg
3271121	Long King Single	1050 mm	2100 mm	200 mm	250 kg
3271125	Long Single	850 mm	2200 mm	200 mm	250 kg
3271126	Long Single	900 mm	2200 mm	200 mm	250 kg
3271128	Long King Single	1050 mm	2200 mm	200 mm	250 kg

Issue Date: 241023
Revision Date: 290726 | V4