

Staebel Opticell AIR1 Reactive Air Mattress

Product Code: 32712

The Staebel Opticell AIR1 reactive hybrid mattress is an advanced hybrid mattress designed to meet the high standards of healthcare and the needs of the end user.



230kg
SWL



High Risk
Mattress

3

Suitable for
Stage 3

T

Treatment
Mattress



Personalised
Pressure Relief

With its foam-free air cells, the Staebel Opticell AIR1 Reactive Hybrid Air Mattress has performance data to support its low pressure, providing excellent pressure relief and comfort (42/44 mm/Hg). The OptiTex™ cover has been tested and complies with being vapour permeability, water resistant, hydrostatic and 4-way stretch.


The OptiCell Air1 is completely silent as it has no electric pump and therefore is well suited for regions with unstable power supply or prone to power outages. The included hand pump is used to check and adjust the air pressure in various sections across the mattress. This mattress combines selected foams to promote immersion and envelopment, key factors in pressure injury prevention. Lastly, the Staebel Opticell AIR1 is fitted with a pressure relief system, meaning there are integrated airlocks that enable the easy disconnection of air cell/s where additional pressure redistribution is required across the surface.

Mattress Features:

- Base foam configuration of high-resilience foam HR50 and hypersoft foam topper. Together, these foams facilitate immersion and envelopment while achieving a high level of comfort for the end user
- The mattress is fitted with a pressure relief system with an integrated airlock that enables the easy disconnection of cell/s
- Staebel OptiCell AIR1 has got documented extra low pressure (42/44 mm/Hg)
- Quality product from Europe
- WARRANTY: 2 Years

Cover Features:

- The OptiTex™ cover can be cleaned with disinfectant and can be machine washed at 95°C. The cells and hoses can be rinsed clean with warm water (max 70°C) or with disinfectant
- The cover surface is made of bi-elastic double-coated textile with soft elastic PU
- 4-way stretch cover material to minimise effects of shear and friction
- The bottom cover has anti-glide protection with four fixing straps to facilitate handling of the mattress
- Drip protection, waterproof zip on four sides. All seams are welded for infection risk mitigation

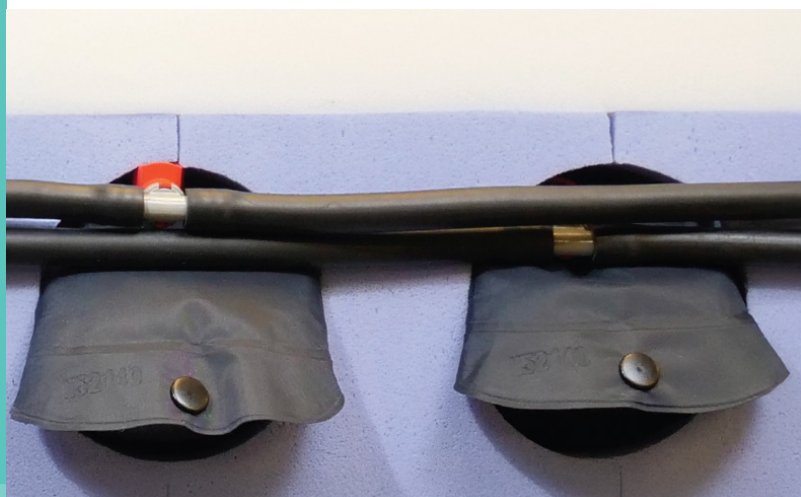

Specifications:

Product Code	3271201	3271203
Size	Single	King Single
Width	850 mm	1050 mm
Length	2000 mm	2000 mm
Thickness	180 mm	180 mm
SWL	230 kg	230 kg

Issue Date: 110325
Revision Date: 150725 | V2