

Staebel Soft EuroFlex Hygiene Mattress

Product Code: 327136SEH

The Staebel Soft EuroFlex Hygiene Mattress has been developed as an aid in the treatment of pressure injuries up to stage 3.



180kg
SWL



Very High Risk
Mattress

3

Aids in the Treatment of up
to & including **Stage 3** PI



Highest Infection
Control




No seams
or zippers

The Staebel Soft EuroFlex Hygiene Mattress has been developed as an aid in the treatment of pressure injuries up to category 3. The mattress is designed for pressure redistribution, specifically for protruding body parts such as shoulders, heels and elbows.

The Staebel Soft EuroFlex Hygiene Mattress is made up of three different density foam cores with comfort and pressure redistribution in mind. The Lentex® laminate surface layer is completely devoid of seams and zippers, which provides increased conditions for good hygiene. Additionally beneficial in healthcare environments, the semi-translucent, white Lentex® laminate allows staff the early detection of strikethrough, minimising healthcare-acquired infections. The addition of an outer mattress cover is possible. From an infection control perspective, the Staebel Soft EuroFlex Hygiene Mattress may be well suited for immunocompromised individuals or those with chronic skin conditions.

Product Features:

- White semi-translucent surface layer also makes it easy for staff to detect contamination on the surface of the mattress, which can be quickly cleaned
- The surface layer is liquid-repellent, and the Active Membrane technology allows the mattress a level of breathability while reducing the risk of liquid ingress
- Lentex® Soft Euroflex has been tested by Swerea IVF according to the measurement method SS 876 00 13
- The foam core consists of three layers with carefully selected foam densities to provide comfort and pressure redistribution
- Can be used without an outer cover. The outer cover may be purchased as an addition
- Safe Working Load of 180kg
- Quality product from Europe
- WARRANTY: Adhesion of Lentex Surface Layer - 5 years, Structural Fatigue of the Core - 3 years

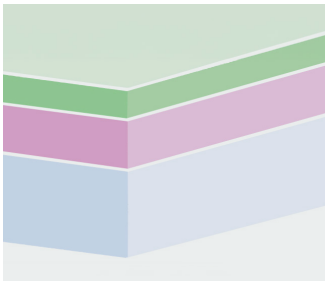

Specifications:

Product Code	3271361	3271362	3271364	3271365
Size	Single	Single	King Single	Euro Wide
Width	850 mm	900 mm	1050 mm	1200 mm
Length	2000 mm	2000 mm	2000 mm	2000 mm
Thickness	130 mm	130 mm	130 mm	130 mm
SWL	180 kg	180 kg	180 kg	180 kg

Issue Date: 121224
Revision Date: 290726 | V4