

CuroCell Area Zone Reactive Air Mattress

Product Code: 36CCAZ

The AREA Zone is a mattress that combines air with foam for effective pressure redistribution and great comfort.



The CuroCell Area Zone Mattress is used as an aid in the prevention and treatment of pressure injuries up to and including category 3. The mattress combines the comfort of the foam with the pressure-relieving properties of the air. The mattress contains nine air cells divided into zones, adapting to the individual and providing a individualised setting. The AREA Zone can be used in all types of care environments and is easy to handle and maintain. The mattress is constructed with multiple layers of high-resilience cold foam. The built-in heel function also reduces pressure on the heels. The mattress maintains a good hygienic standard and is manufactured with welded seams.



230kg
SWL




High Risk
Mattress

3

Suitable for
Category 3


P

Prevention
Mattress

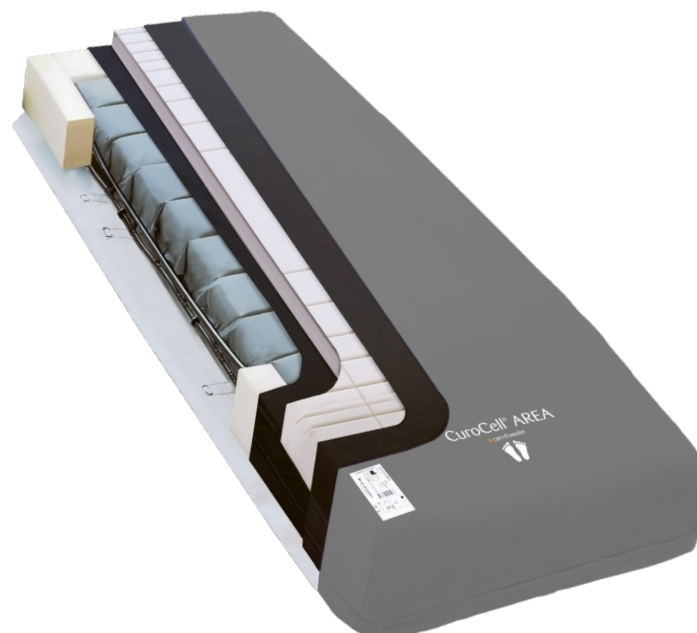

Product Features:

- Automatic control unit which adapts individual to the user
- User friendly interface
- Audible and visual notifications
- Three program modes: pulsating mode, alternating mode and constant low pressure
- Silent running and vibration free
- Air flow control; Constant temperature, good user comfort, energy reduction
- Soft foam surface
- Smart cable management
- Easy handling and cleaning
- Heel slope, for pressure redistribution of heels
- CPR-function with transport function

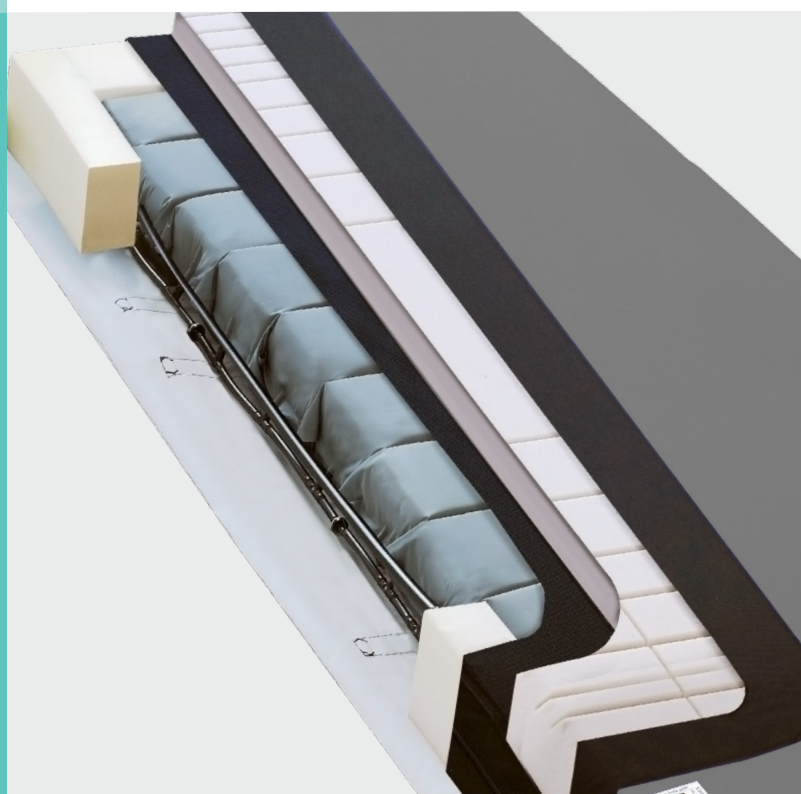

Specifications:

Product Code	3678Z	3679Z	3686Z	3680Z
Width	850 mm	900 mm	1000 mm	1050 mm
Length	2000 mm	2000 mm	2000 mm	2000 mm
Thickness	180 mm	180 mm	180 mm	180 mm
SWL	230 kg	230 kg	230 kg	230 kg