

CuroCell A4 Cirrus Active Air Mattress

Product Code: 36CCA4C

An automatic mattress replacement system with foam and air combined.



CuroCell® A4 Cirrus is used as an aid in the prevention and treatment of pressure injuries up to and including category 4 and unclassified pressure injuries⁽¹⁾. It is a replacement system suitable for all healthcare environments.

CuroCell® A4 Cirrus combines air with a soft foam surface and provides good user comfort and pressure redistribution⁽³⁾. The added heel function offers pressure redistribution on the heels.

CuroCell® Cirrus has an automatic control unit with a user-friendly interface. The mattress system regulates independently and without manual adjustment the inner pressure of the mattress to different values based on the user's weight and position. No manual action needs to be done to adjust the internal pressure of the mattress.



230kg
SWL



High Risk
Mattress

4


Suitable for
Category 4

P

Prevention
Mattress

Product Features:

- Automatic control unit which adapts individual to the user
- User friendly interface
- Audible and visual notifications
- Three program modes: pulsating mode, alternating mode and constant low pressure
- Silent running and vibration free
- Air flow control; Constant temperature, good user comfort, energy reduction
- Soft foam surface
- Smart cable management
- Easy handling and cleaning
- Heel slope, for pressure redistribution of heels
- CPR-function with transport function

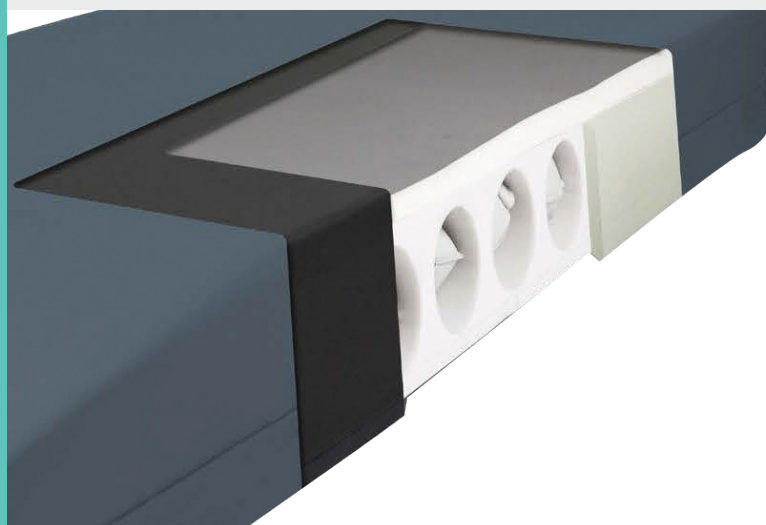

Specifications:

Product Code	3674A4	3675A4	3687A4	3676A4
Width	850 mm	900 mm	1000 mm	1050 mm
Length	2000 mm	2000 mm	2000 mm	2000 mm
Thickness	210 mm	210 mm	210 mm	210 mm
SWL	230 kg	230 kg	230 kg	230 kg