

CuroCell® Area Zone Reactive Hybrid Air Mattress

Product Code: 36CCAZ

The AREA Zone is a hybrid mattress that combines air with foam for effective pressure redistribution and great comfort.



230kg
SWL



High Risk
Mattress

3

Aids in the Treatment of up
to & including **Stage 3** PI




Integrated Heel
Function/Slope



4-Way
Stretch Cover


The CuroCell® Area Zone Reactive Hybrid Mattress is used as an aid in the prevention and treatment of pressure injuries up to and including stage 3. The mattress combines the comfort of the foam with the pressure-relieving properties of air, containing 9 air cells divided into zones, adapting to the individual. The valve system allows air flow between the chambers, which in turn provides a unique, custom low-pressure profile for any patient. Upon ingress, a natural chain of events begins in each of the connected air chambers. The slow-release valves allow for pressure changes to take place very gently.

The CuroCell® Area Zone is designed with an inner anti-shear layer between the outer cover and inner foam, whereby shear is one of the leading known causes for pressure injury development. The CuroCell® Area Zone can be used in all types of care environments and is easy to handle and maintain, making it a suitable option if looking to reduce mattress maintenance or in regions where the power supply is unstable and of concern.

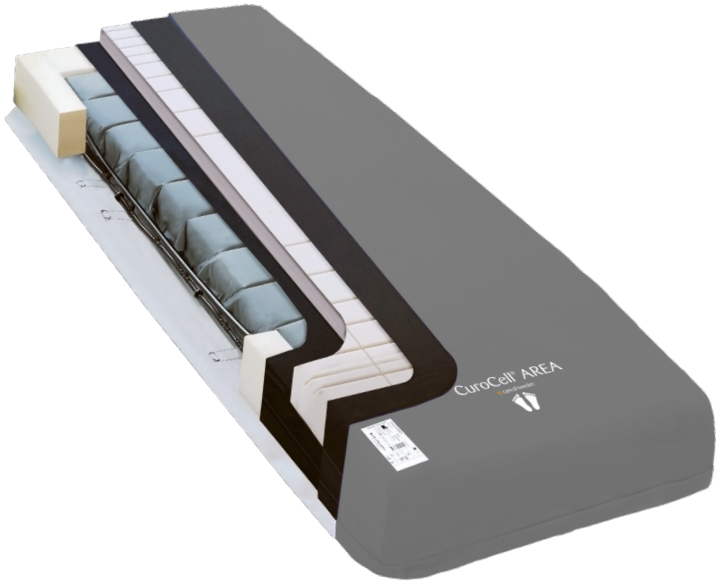

Product Features:

- Suitable for high-risk individuals
- Static air replacement mattress system (not an overlay)
- Reinforced sides for improved lateral stability and to assist with independent bed transfers
- Unique integrated heel protection area/heel slope
- Anti-shear inner cover
- 4-way stretch cover material to minimise the effects of shear and friction
- Cover Material is made from 61% polyamide, 39% polyurethane coating
- No pump or power cord for minimal maintenance
- Safe Working Load 230kg
- Quality product from Europe
- WARRANTY: 5 years

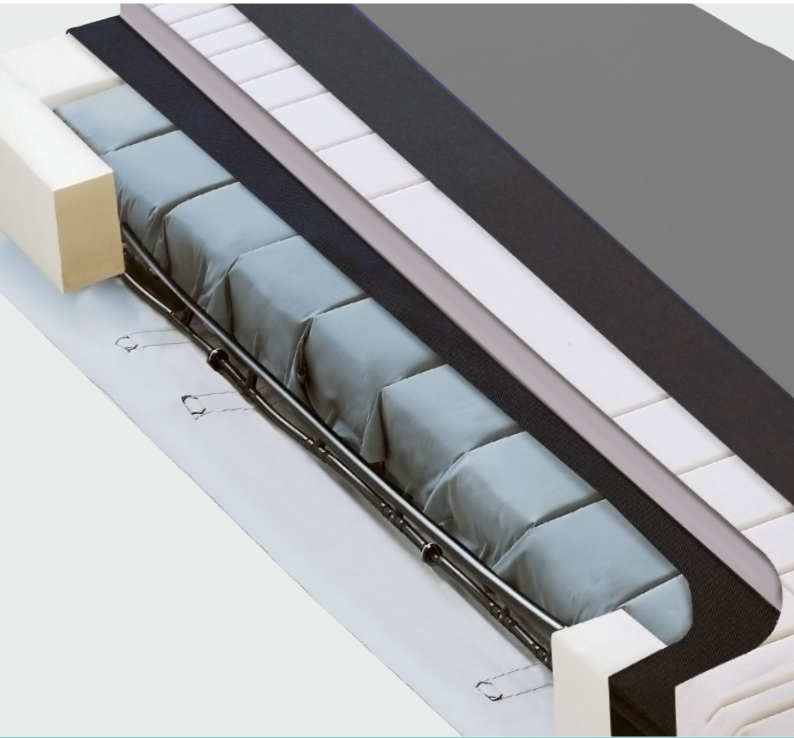

Specifications:

Product Code	3678Z	3679Z	3686Z	3680Z	3681Z
Size	Single	Single	King Single	King Single	Extra Wide
Width	850 mm	900 mm	1000 mm	1050 mm	1200 mm
Length	2000 mm	2000 mm	2000 mm	2000 mm	2000 mm
Thickness	180 mm	180 mm	180 mm	180 mm	180 mm
SWL	230 kg	230 kg	230 kg	230 kg	230 kg

Issue Date: 150725
Revision Date: 290726 | V3