

Staebel Eurocare Foam Mattress

Product Code: 3728

Promoting pressure reduction, assisting in the prevention of pressure injuries while enhancing patient well-being.



25 - 150kg
SWL



Low Risk
Mattress



Prevention
Mattress



Firm
Surrounds



Multi-Way
Stretch Cover


The Staebel EuroCare Foam Mattress is a premium-quality solution designed to deliver exceptional support and comfort for patients in home-care and nursing environments. Engineered with a profile-cut pressure-relieving foam insert, this mattress promotes effective pressure reduction, assisting in the prevention and management of pressure injuries. This mattress is considered a comfortable yet cost-efficient solution for both home care and residential aged care settings and is a reversible support surface suitable for low-pressure injury risk individuals.



Product Features:

- **Castellated foam surface:** to enhance envelopment and provide pressure redistribution
- **Firm surrounds:** enabling independent bed transfers and mobility
- **Profile-Cut Foam Insert:** Strategically designed to reduce pressure points, ensuring optimal weight distribution for patients at risk of pressure injuries.
- **Enhanced Comfort:** Provides a comfortable sleeping surface that supports both patient health and restfulness.
- **Durability and Practicality**
- **Reinforced Edges:** Constructed with unique edge reinforcements to enhance durability and prolong mattress lifespan.
- **Waterproof Zip:** Helps maintain a clean and sanitary environment by providing a barrier against fluid ingress.
- **Vapour-Permeable, Water-Resistant Cover:** The mattress cover balances breathability with protection, ensuring comfort while preventing moisture buildup.
- **Fire-Retardant Material:** Meets stringent safety standards for fire resistance, offering peace of mind in care settings.
- **Welded Seams:** Precision-sealed seams enhance infection control by reducing potential contamination points.
- **Inner Anti-Shear Cover:** Reduces friction and shear, which are key contributors to pressure injury development
- **WARRANTY:** Foam core - 5 years , Cover - 2 years

Available in custom sizes - Contact the team now

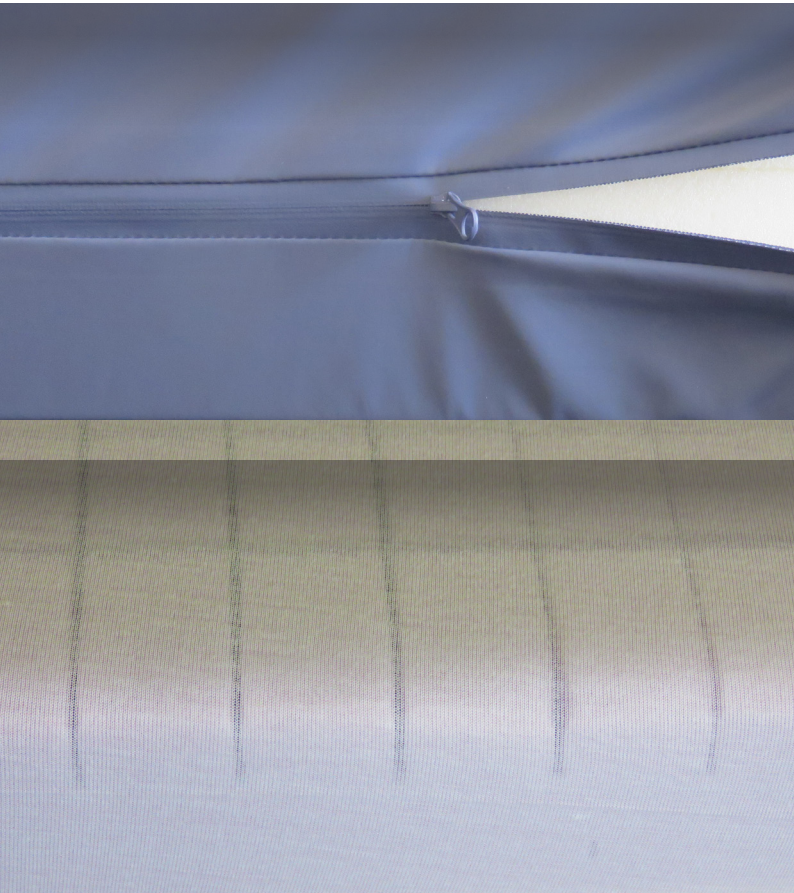

Specifications:

Product Code	3728	3729	3728KS	3728D	3728Q	3728K
Size	Single	Wide	King Single	Double	Queen	King Double
Width	850 mm	950 mm	1050 mm	1350 mm	1500 mm	1800 mm
Length	1980 mm	1980 mm	1980 mm	1980 mm	2000 mm	2000 mm
Thickness	150 mm	150 mm	150 mm	150 mm	150 mm	150 mm
SWL	25 - 150 kg	25 - 150 kg	25 - 150 kg	25 - 150 kg	25 - 150 kg	25 - 150 kg

Issue Date: 161224
Revision Date: 130525 | V3