

# Staebel Zonecare 300 Bariatric Medical Mattress

Product Code: 3M9300

The Staebel Zonecare 300 Bariatric Medical Mattress is a cost-efficient pressure care solution, purpose-built to support bariatric users with enhanced comfort, stability, and safety.



300kg  
SWL



Low Risk  
Mattress



Prevention  
Mattress



Bariatric  
Use



Firm  
Surrounds

Engineered with a triple-layered foam core, the Staebel Zonecare 300 Bariatric Medical Mattress delivers a combination of pressure redistribution, postural support, and durability. A visco-grade upper for immersion, a mid-core for comfort and support, and a firm base to prevent bottoming out during transfers and repositioning.

The Staebel Zonecare 300 Bariatric Medical Mattress has an anti-microbial treatment applied to the foam, enhancing hygiene and core longevity, while strengthened sides increase structural integrity. The mattress is encased in a durable, vapour-permeable multiway stretch cover, reducing shear forces and maintaining skin integrity. Fire retardant and compliant with key safety standards, the Staebel Zonecare 300 Mattress is a reliable choice for bariatric care in demanding clinical environments.

## Product Features:

### VISCO-GRADE, UPPER IMMERSION LAYER :

- Alleviates localised pressure
- Provides maximum contact with the patient for the distribution of interface pressure
- Visco-elastic cell structure provides a low-shear surface and immersion

### TRIPLE LAYERED CORE:


- Premium Grade Mid Level, provides open-cell high quality foam
- Provides exceptional comfort and postural support for repositioning and patient care
- Foam conforms to AS/NZS 4088.1 (Int):1993 and FMVSS302 & SAE T2002
- Anti-microbial treatment for maximum longevity

### FIRM BASE LAYER:

- Protects against 'bottoming-out' during transfer, repositioning and whilst in the inclined/ fowler position
- Provides exceptional comfort, support for repositioning and patient care
- Premium quality, high-performance foam

### Cover Features:

- Multiway stretch for shear protection
- Vapour permeable and breathable
- Flexible and gentle to the skin
- Fire retardant, meets BS 6807 ignition source 5 and BS 7175 ignition source 5 standards
- WARRANTY: Foam core - 5 years , Cover - 1 years

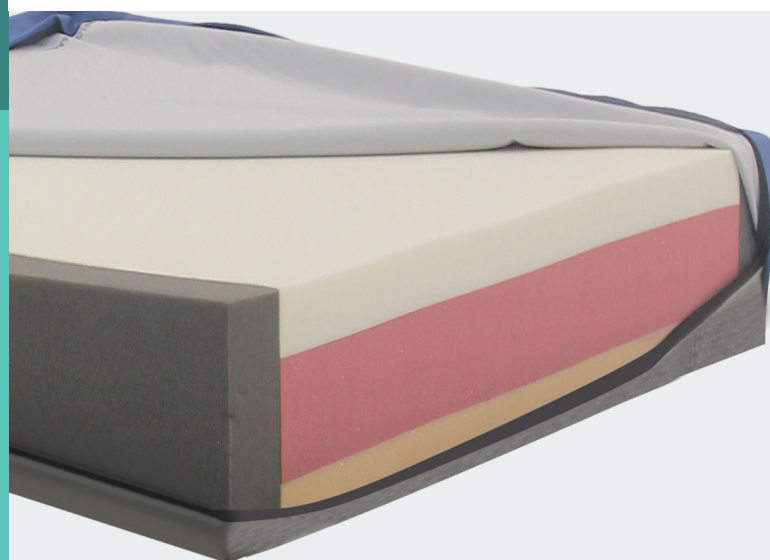

## Specifications:

Product Code	3M9300-KS	3M9300-EW	3M9300-D	3M9300-Q	3M9300-K
Size	King Single	Extra Wide	Double	Queen Double	King Double
Width	1050 mm	1150 mm	1350 mm	1480 mm	1780 mm
Length	1980 mm	1980 mm	1980 mm	1980 mm	1980 mm
Thickness	150 mm	150 mm	150 mm	150 mm	150 mm
SWL	300 kg	300 kg	300 kg	300 kg	300 kg

Issue Date: 231224  
Revision Date: 190326 | V2